

Adenoviral Knockdown Virus Data Sheet

Vector Information

Vector: E1/E3 deleted adenoviral serotype 5 vector
 Expression cassette: shRNA under control of human U6 promoter

Target Information

Target Gene: bcl2
 Specimen: Homo sapiens
 RefSeq#: NM_000657

Validation

Validation Method: qRT-PCR, Relative Quantification
 Cell line used for validation: SK MEL 28
 RNA prepared: 48h post transduction
 Total RNA used in RT: 1 µg
 cDNA Primer: Random Hexamer
 Real Time Platform: Light Cycler 480 (Roche)
 Detection Method: SYBR I
 Reference Gene: ACTB
 Amplicon Integrity verified by: Melting Point Analysis

Results

Viral vector	Target mRNA expression rel. To Ctrl
Ad 095	38%
Ad 126	20%
Ad 005	27%
Ad 094	28%

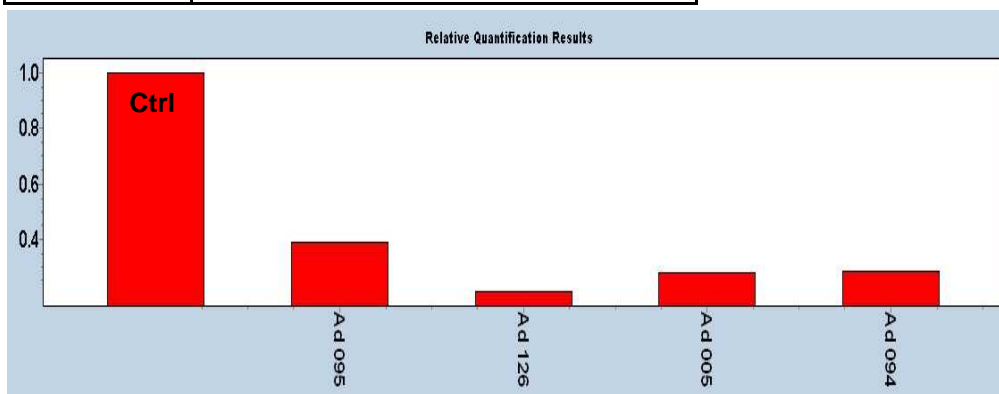


Figure 1

Reduction of bcl2-mRNA 48h after transduction with adenoviral vectors expressing 4 different target-specific shRNAs in comparison to cells transduced with control vector (Ctrl)